

ACP-100 MIRROR CHART

SPECIFICATIONS

Charts.....	32 Visual targets
Selection Speed.....	Average 0.45 sec
Projection distance.....	1.1m
Auto-off function.....	5 minutes
Lamp.....	12V/35W
Power source.....	Main body: AC110V/AC220V 50/60Hz Remote control:DC9V
Power consumption.....	50VA
Dimensions.....	Main body: 400X280X560mm Remote control:248X70X18mm
Net weight.....	Main body: 20kgs Remote control: 0.22kg



NV-100 NEAR VISION TESTER

FEATURES

The refractive tests include:

- Measurement of visual acuity using Snellen letter and illiterate E-chart,and Times Roman reading types.
- Accurate determination of presbyopic additions and other near corrections.
- Near vision astigmatism tests.
- Amplitude of accommodation including monocular measurements under binocular conditions.

The Binocular balance tests include:

- Fixation disparity test for uncompensated heterophoria motor eye dominance. Accommodation-convergence relationship.
- Quantitative evaluation of foveal suppression in heterophoria patients.
- Near point of convergence.
- Targets for subjective cover tests for heterophoria and harmonious ARC.

SPECIFICATIONS

Test Distance.....	36cm
Distant Control.....	60cm
Size.....	220(H)×210(W)×42(D)mm
Weight.....	780g
Power Supply.....	LR6 AA/1.5V×3

The NV-100 near vision tester enable a complete refractive and binocular vision investigation to be carried out with the head and eyes in their normal posture at any working distance. It's in English. But French and Dutch version are also available.



NEAR VISUAL ACUITY CHART

FEATURES

- Double side display(2 charts). Amsler grid is the standard chart for the other side.
- 1100mAh rechargeable Li-ion battery offer longer working hour.
- The angle of the base can be adjusted.


SPECIFICATIONS

Test distance.....	30cm to 40cm
Dimension.....	196X145X10mm
Power supply.....	AC adaptor, rechargeable Li-ion battery

STANDARD ACCESSORIES:




More charts of near vision chart

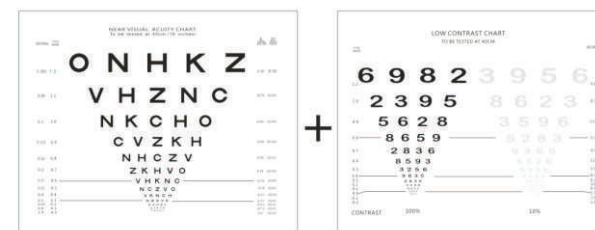

Letter chart

Low contrast

NV-194

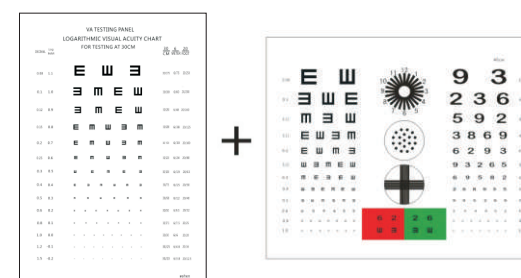

E chart

E and function chart

NV-191




Vertical chart

Lateral chart

NV-195

NV-145

AMSLER GRID
back side of NV195
and NV 145

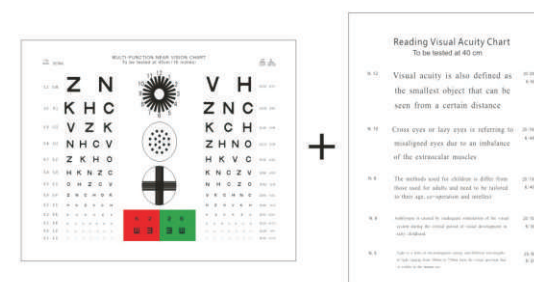

Letter and function chart

Reading chart

NV-134

